



50/10 GENESI M

Cod. D-GE150-040-M


Oil free compressor, FAD at 5 bar 135 l/min, tank capacity 40 L, with membrane air dryer and 0,01 micron filtration.

- Maximum simultaneous use: 1
- 🔄 Dental Chairs: 2

Medical Device CE 0425

			ē	5		(٧	Ļ,	\bigcirc		kg
Air Dryer	Air capacity	FAD Fad @ 5bar	Motor	Voltage Hz	Receiver	Noise Level	Duty Cycle	max P.	Size	net Weight
Yes	m³/min l/min 0.182 182	m³/min l/min 0.135 135	kW 1.1	230V 50Hz	l 40	db(A) 65	50 %	bar 7	mm 710x410x770	Kg 42

Air delivery performance

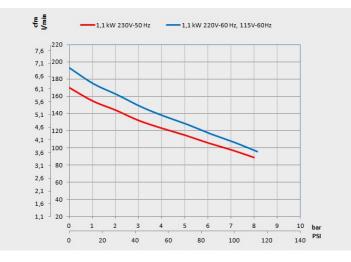



Maintenance kits

Kit.	Description	Туре	Hours
TA-FL-024-K-SCC	GENESI M ordinary maintenance kit Ordinary maintenance kit for GENESI compressors with dryer. The kit consists of the cartridges of the dryer filtration battery (2 cartridges: 5 + 0.01 micron) and the relative O-rings, the pumping element suction filter and the non-return valve pad. Annual replacement of kit components is recommended. Maintenance operations can be performed at a higher frequency if the environment in which the compressor works is particularly dirty.	Ordinary	1.000h/ly
MG-GC-011-K	Extra ordinary maintenance kit for GENESI S and M versions The GENESI extraordinary maintenance spare parts kit includes: a suction filter, two elastic rings, O-rings, a line filter cartridge and two complete valve plates. The recommended maintenance intervals for this compressor are 1500 hours / 2 years. Maintenance operations can be performed at a higher frequency if the environment in which the compressor works is particularly dirty.	Extraordinary	1.500h/2y